



quality dental supplies



gifting smiles worldover

auxillaries



Etching Gel

DESCRIPTION AND FEATURES

- 37% phosphoric acid gel
- Silica free, high quality, non - drying, medium viscosity gel
- Thixotropic ,stays put, does not run
- Visible blue colour
- Has excellent water solubility
- Disposable tips save time and are more sanitary

INDICATIONS FOR USE

- Before use of bonding material
- Before use of filling materials
- Before use of orthodontic brackets

TECHNICAL DATA

Excellent water solubility	
pH	2
Viscosity	Thixotropic
Shelf life	2 years



PACKAGING

- 1 syringe of 3gm etch
- 2 Applicator tips

*Also available in 10 ml dropper bottle pack

EDTA Gel

DESCRIPTION AND FEATURES

- Acts as a chelating agent for dentin calcium
- Lubricates the root canal and facilitates endodontic instrument penetration
- Makes shaping more efficient
- Superior wettability ensures that it gets uniformly coated all over the canal and can be washed off easily

INDICATIONS FOR USE

- Facilitates root canal preparation

PACKAGING

- 1 Syringe of 3gm EDTA Gel
- 2 Applicator tips



- etching gel
- RCEZ gel
- EZ white
- G.P. Solve
- zinc oxide
- eugenol
- varnish
- conditioner
- cal - plus

www.dtechasia.com



quality dental supplies



gifting smiles worldover

auxillaries



EZ White Gel

DESCRIPTION AND FEATURES

- The viscosity of the gel provides maximum contact with tooth surface
- Quick & effective whitening
- EZ-White will remove most of the commonly found stains quickly & effectively
- Full treatment lasts just few days
- Fresh Mint flavour
- Available in 11%, 16% & 21% concentrations



INDICATIONS FOR USE

- For bleaching/whitening of teeth


PACKAGING

- 1 Syringe of 2.5gm
- 2 Applicator tips

Cal - plus

DESCRIPTION AND FEATURES

- Radio opaque Calcium hydroxide paste
- Ready to use, water based paste
- Highly alkaline

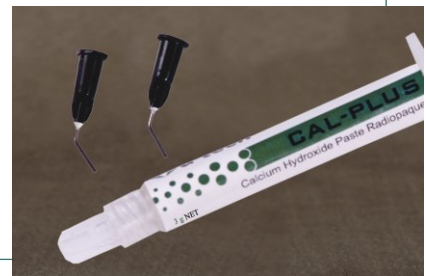

INDICATIONS FOR USE

- Pulp capping
- As a cavity base
- As a neutralizing agent or protecting dentine against chemical or thermal reaction

PACKAGING

- Syringe of 3gm/ 2 Applicator tips



- etching gel
- RCEZ gel
- EZ white
- G.P. Solve
- zinc oxide
- eugenol
- varnish
- conditioner
- cal - plus

Eugenol

DESCRIPTION AND FEATURES

- A brownish viscous liquid with concentrated eugenols
- Best grade used



PACKAGING

- 1 Bottle of 10 ml/ 110ml

Zinc Oxide

DESCRIPTION AND FEATURES

- A chemically pure, arsenic free material popular for its fine quality
- Best grade used

INDICATIONS

- As a temporary filling material to be used in conjunction with Eugenol.

PACKAGING

- 1 Bottle of 110 gms



www.dtechasia.com



quality dental supplies



gifting smiles worldover

auxillaries



GP-Solv

DESCRIPTION AND FEATURES

- Does not contain chloroform
- Refreshing Citrus Fragrance
- Contains natural non-hazardous solvent

INDICATIONS FOR USE

- Root Canal Filling Remover
- Effectively softens GP-Points & Zinc Oxide eugenol cement



PACKAGING

- 1 Bottle of 10 ml

Varnish

DESCRIPTION AND FEATURES

- For Protective Coating on Glass Ionomer restoration
- Protects from moisture during the initial setting
- Does not contain acetone or chloroform
- Contains natural non-hazardous solvent
- Refreshing Citrus fragrance

INDICATIONS FOR USE

- Protection of Glass Ionomer restorations from moisture during initial setting period.



PACKAGING

- 1 Bottle of 10 ml

Conditioner

DESCRIPTION AND FEATURES

- D-tech Dentine conditioner is a mild solution of polyacrylic acid
- Removes the smear layer from dentin and improves the bonding of Glass Ionomer restorations and linings

INDICATIONS FOR USE

- The D-Tech Dentine conditioner effectively removes the smear layer & significantly improves the bonding of Glass Ionomer to dentin.



PACKAGING

- 1 Bottle of 10 ml



- etching gel
- RCEZ gel
- EZ white
- G.P. Solve
- zinc oxide
- eugenol
- varnish
- conditioner
- cal - plus

www.dtechasia.com